





Evolve Cemented Stem



Evolve cemented stem



Launched in 2010

Over 4,000 supplied worldwide

FDA approved

CE marked

3A* ODEP rated







Evolve cemented stem




Cemented Hip System with proven geometry and material - High Nitrogen Stainless Steel ISO5832-9







Evolve features

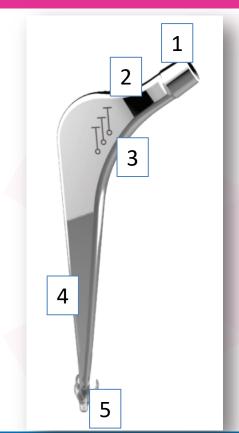

1. 12/14 taper

2. Multiple head offsets

3. Collarless neck facilitates intraoperative adjustment

 4. Highly polished double-taper helps create compressive, radial loading.
Reduces friction between the stem and cement mantle

5. PMMA distal centraliser reduces pointloading on the distal cement mantle








Sizing table

Evolve Masters					
P/N	Size	LT	0	a°	G
111-242-351	35.5 Size 1	127	35.5	125	10.9
111-242-370	37.5 Size 0	150	37.5	125	10.4
111-242-371	37.5 Size 1	150	37.5	125	10.9
111-242-372	37.5 Size 2	150	37.5	125	11.4
111-242-373	37.5 Size 3	150	37.5	125	11.9
111-242-440	44 Size 0	150	44	125	10.4
111-242-441	44 Size 1	150	44	125	11.0
111-242-442	44 Size 2	150	44	125	11.4
111-242-443	44 Size 3	150	44	125	11.9
111-242-444	44 Size 4	150	44	125	12.6
111-242-501	50 Size 1	150	50	125	10.9
111-242-502	50 Size 2	150	50	125	11.4
111-242-503	50 Size 3	150	50	125	11.9
111-242-504	50 Size 4	150	50	125	12.6

Evolve					
111-162-133	SO Size 1	116	37.5	130	11.0
111-162-142	HO Size 1	116	40.5	125	11.0
111-162-240	SO Size 2	125	39.5	130	11.0
111-162-244	HO Size 2	125	44	125	11.0
111-162-342	SO Size 3	134	41.5	130	11.0
111-162-364	HO Size 3	134	46	125	11.0
111-162-444	SO Size 4	143	43.5	130	11.5
111-162-448	HO Size 4	143	48	125	11.5
111-162-546	SO Size 5	151	45.5	130	12.0
111-162-550	HO Size 5	151	50	125	12.0

		Ø28	Ø32	Ø36	Ø40
s	Croanic Croanic	-15 (11-152-01) (11-152-01)	-4.0 (111-152-62) (111-152-62)	-4.0 (111-152-411 (111-152-411	-4.0 (111-152-64) (111-152-64)
м	OCo Cetanik	0.00 111-152-012 111-152-612	.0 111-152-022 111-152-022	0.00 111-152-012 1111-152-012	.0 111-152-002 111-152-002
L	CiCo Cesanic	+4.0 (111-152-013 (111-152-013	+4.0 111-152-623	+4.0 (111-152-413) (111-152-413	+4.0 (11-152-04)
хL	Crownia		+8.0 111-152424	+8.0 111-152-634 111-152-634	+8.0 (111-152-644

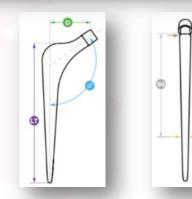








6

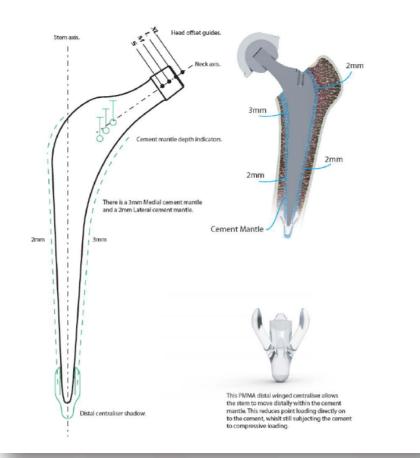










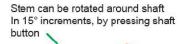


Soft release: Pulling trigger removes locking pin and allows inserter to rotate off trunnion collar without disturbing the stem







Instrument tray

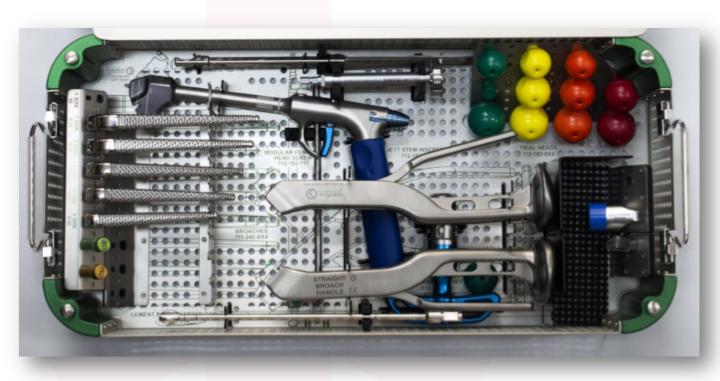








Instruments

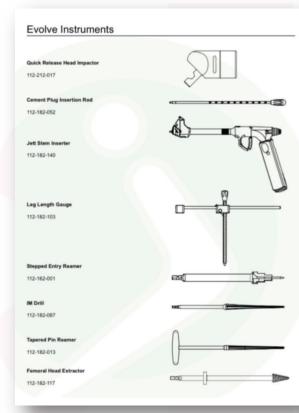

Straight Broach Handle

112-182-112

Straight Broach Handle

112-182-050

Short Broach Handle 112-25-0146

Dual Offset Broach Handle



Curved Broach Handle

112-182-001

Curved Broach Handle

112-182-166


Calcar Reamen

P12-01-0012

Osteotomy Guide

112-212-005

Quick Release Box Chisel

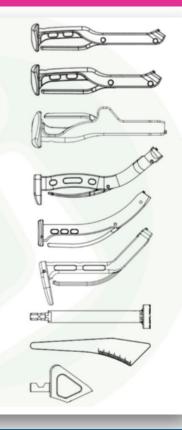








Instruments

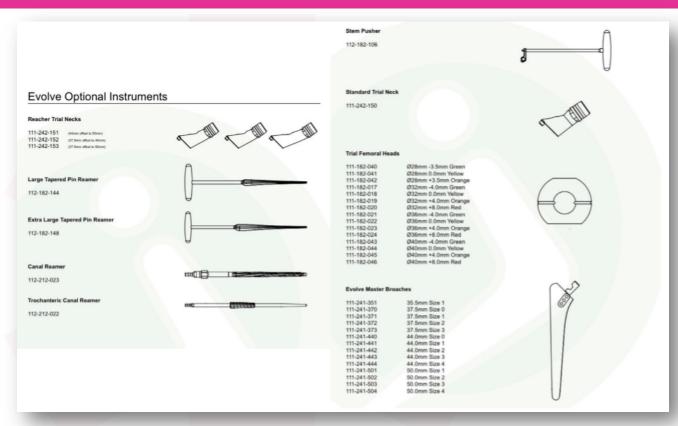








X-ray templates

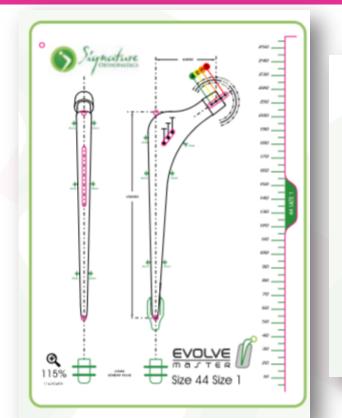

Evolve Templates

111-242-300 111-242-302 111-242-303 111-242-303 111-242-305 111-242-306 111-242-306 111-242-307 111-242-308 111-242-309 111-242-310 111-242-311 111-242-312 111-242-313 35.5mm Size 1 37.5mm Size 0 37.5mm Size 1 35.5mm Size 2 37.5mm Size 2 37.5mm Size 3 44.0mm Size 0 44.0mm Size 1 44.0mm Size 1 44.0mm Size 3 44.4mm Size 4 50.0mm Size 1 50.0mm Size 2 50.0mm Size 3







Evolve implants

Evolve Master Implants

111-242-351	35.5mm Size 1
111-242-370	37.5mm Size 0
111-242-371	37.5mm Size 1
111-242-372	35.5mm Size 2
111-242-373	37.5mm Size 3
111-242-440	44.0mm Size 0
111-242-441	44.0mm Size 1
111-242-442	44.0mm Size 2
111-242-443	44.0mm Size 3
111-242-444	44.4mm Size 4
111-242-501	50.0mm Size 1
111-242-502	50.0mm Size 2
111-242-503	50.0mm Size 3
111-242-504	50.0mm Size 4
Evolve Stems	

Distal Centraliser	
111-182-002	Wide
111-162-003	Narrow
Cement Plug	
112-182-150	14mm
112-182-151	17mm
112-182-152	21mm
112-182-153	25mm
112-182-150 112-182-151 112-182-152	17mm 21mm

11-162-133	Size 1 Standard Offset
11-182-142	Size 1 High Offset
11-162-240	Size 2 Standard Offset
11-182-244	Size 2 High Offset
11-162-342	Size 3 Standard Offset
11-162-364	Size 3 High Offset
11-182-444	Size 4 Standard Offset
11-162-448	Size 4 High Offset
11-182-548	Size 5 Standard Offset
11-182-550	Size 5 High Offset









Thank you

